

# 苏州市相城区消毒站

## Suzhou City Xiangcheng District Sterilization Station

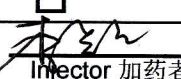


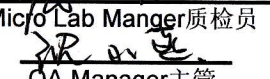
### 环氧乙烷灭菌批记录(EO-Sterilization Batch Report)

Furnace No箱号: Y1-01

Sterilization No.灭菌批号: 0601

Date 日期: 2019-06-28

Article Name 物品名称	Production No 生产批号	Expiration Date 有效期	Specification 规格/型号	Quantity 数量	Remark 备注	Parameters During Sterilization Cycle 灭菌过程参数	Required 设定值	Actual 实际值
Cartridge Needle	0601	06/2024	1RL 0.18MM、 3RL 0.22MM、 3RS 0.35MM etc	29CTNS	CERUS	灭菌温度Chamber Temperature while Sterilization(°C)	47-53	S:48.5 E: 52.6
						灭菌压力Chamber Pressure while Sterilization(kPa)	-30-20	S:-27.1 E:-15.5
						灭菌湿度Chamber Humidity while Sterilization(%RH)	40-80	33.5-52.8
						灭菌时间Exposure time(min)	480	480
						预真空度Pre-Vacuum Value(kPa)		
						加湿压力幅度Pressure Amplitude while Humidifying(kPa)		
						加湿前后压力Pressure Changed while Humidifying(kPa)		
						加药压力幅度Pressure Amplitude while EO Injection(kPa)		
						加药前后压力Pressure Changed while EO Injection(kPa)	//	S:-69.2 E: -35.2
						加药重量Weight of EO used(kg)	3.8-4.2	12

Sterilization Phase 灭菌过程	Start Time 开始时间	End Time 结束时间	Actual Time 实际时间	Sterilization Information 灭菌情况	Release 验收
预热 Conditioning	9: 17	13: 17	240	Normal 正常	<input checked="" type="checkbox"/> <b>Passed</b> 合格 <input type="checkbox"/> <b>Rejected</b> 拒收 <input type="checkbox"/> <b>Hold</b> 待定  Inspector 加药者  EO Manger 灭菌主管  Micro Lab Manger 质检员  QA Manager 主管
抽真空 Air Removal	13: 17	13: 24	7	Normal 正常	
保压时间 Pressure Time	13: 24	13: 29	5	Normal 正常	
充入EO EO Injection	13: 35	13: 50	15	Normal 正常	
灭菌时间 Exposure Time	13: 50	21: 50	480	Normal 正常	
换气时间 EO Removal	21: 50	22: 16	26		
清理完毕 Finishing	22: 16	22: 48	32		

